


Additional Tissue Expanders



sientra®

Calf Tissue Expander


Smooth Surface
Remote Port

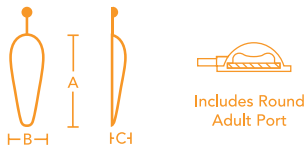

DIMENSIONS				REFERENCE
A cm	B cm	C cm	VOL cc	
19.2	5.3	2.6	145 cc	880-145
22.9	5.7	2.6	185 cc	880-185


Rectangular Tissue Expander

Smooth Surface
Remote Port

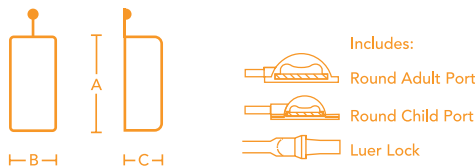

DIMENSIONS				REFERENCE
A cm	B cm	C cm	VOL cc	
5.1	2.5	2.5	30 cc	810-030
3.8	3.8	3.8	50 cc	810-050
6.0	4.6	5.3	100 cc	810-100
8.7	4.1	3.3	105 cc	810-105
11.0	5.1	3.8	160 cc	810-160
6.0	5.0	7.0	200 cc	810-200
9.5	5.5	4.0	205 cc	810-205
9.6	5.9	5.0	250 cc	810-250
14.1	6.2	4.7	320 cc	810-320
13.6	5.5	5.6	400 cc	810-400
15.8	6.9	5.1	500 cc	810-500
18.0	5.0	7.0	600 cc	810-600
18.0	8.0	5.0	640 cc	810-640
28.0	5.0	5.0	700 cc	810-700
15.2	7.3	6.3	720 cc	810-720
18.5	8.7	5.7	800 cc	810-800
20.2	8.0	7.0	1000 cc	810-1000

Smooth Surface
Integrated Port

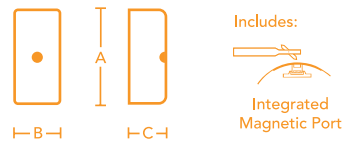

DIMENSIONS				REFERENCE
A cm	B cm	C cm	VOL cc	
11.0	5.1	3.8	160 cc	819-160
6.0	5.0	7.0	200 cc	819-200
9.5	5.5	4.0	205 cc	819-205
9.6	5.9	5.0	250 cc	819-250
14.1	6.2	4.7	320 cc	819-320
13.6	5.5	5.6	400 cc	819-400
15.8	6.9	5.1	500 cc	819-500
18.0	5.0	7.0	600 cc	819-600
18.0	8.0	5.0	640 cc	819-640
28.0	5.0	5.0	700 cc	819-700
15.2	7.3	6.3	720 cc	819-720
18.5	8.7	5.7	800 cc	819-800
20.2	8.0	7.0	1000 cc	819-1000



Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.


Tissue Expander Ports

The use of a 21-gauge needle or smaller is recommended.

PORT OPTIONS:

 = Round Adult Port

 = Round Child Port

 = Luer Lock

 = Integrated Magnetic Port



ROUND ADULT PORT:

Base Diameter 34 mm



ROUND CHILD PORT:


Base Diameter 22 mm

The base of each Sientra Round Port contains a puncture-resistant, radiopaque, stainless steel disk. The disk prevents the needle from penetrating the base when the needle is used as directed (inserted perpendicular to the port and pushed all the way to the bottom).




LUER LOCK:

This luer accompanies intraoperative expanders and some expanders for slow expansion as an external port.

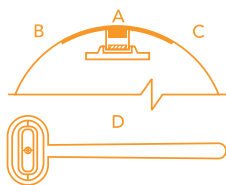


INTEGRATED MAGNETIC PORT FOR SLOW EXPANSION WITH A MAGNETIC FINDER:

This port has a Silicone Shell Septum (A), reinforcement to prevent folds (B), magnetic disk with stabilizer (C), and comes with a magnetic finder (D).

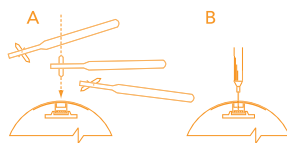

LOCATING THE PORT:

(A) Mark the exact place for the injection.
(B) When introducing the needle, make sure it passes through the septum completely

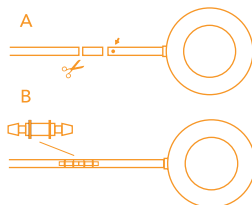

ASSEMBLY SEQUENCE:

(A) Cut the tube to the required length without going past the mark. (B) Place expander tube into the free end of the connective pin until it meets the other tube.

IMPORTANT: Whenever necessary, we can supply any kind of expander with a Child or Adult Round Port or Luer Lock.



Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.

sientra.

PO Box 1490 | Santa Barbara, CA 93116-1490 | Tel: 805.562.3500 | Toll Free: 888.708.0808
Fax: 805.562.8401 | www.sientra.com

MDC-0026 R2